


JMA

Disposable tips



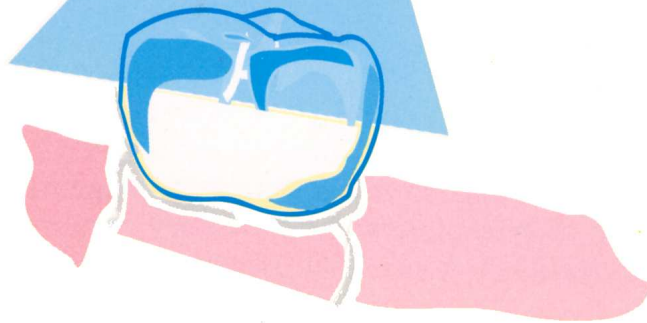


In the fight against the risks of the crossing contamination in a dental office, the cleaning and the disinfection of the multifunctional syringe remain a true problem.

Although the syringe can be sterilized with the autoclave, time does not enable you to sterilize it between each patient.


Concept JMA, thanks to an adapter autoclavable and to a disposable tip, brings the solution at your problem.

JMA Concept is formed of two components:

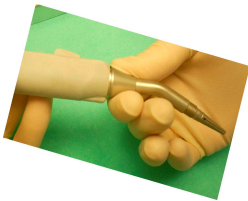

- **The adapter** adapts manually to the original multifunctional syringe. There are several references of the adapters, depending on a syringe used by the practitioner.

- **The cannula**, disposable tip fit together to the end of the adapter by hand screwing (3/4 of turn).

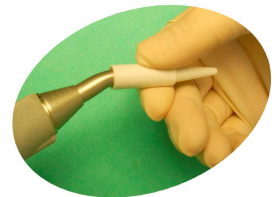

Patented system

JMA Protection

- ✓ The cannula is for a single use
- ✓ The adapters can be autoclaved



- ✓ The objective of JMA system to eliminate any risk of contamination for you and your patient
- ✓ The cannula must be discarded and replaced after each patient

- ✓ The tips are packed in a heat-welded plastic sachet



- ✓ The plastic sachet is conditioned in a carton sealed by a label

- ✓ The adapters are protected by a fine foam plate and are conditioned in a rigid plastic box



JMA Money and Time Saving

- ✓ JMA System represents a considerable time saving for the cleaning and the disinfection of your syringe

- ✓ Conditioning by box of 400 non sterile tips available in white or green



- ✓ JMA cannulas are made from entirely recyclable food plastic
- ✓ JMA adapters are guaranteed for one year

JMA Performance

✓ In a straight line water jet allowing to reach any point in mouth precisely


✓ Slightly bent adapter allowing a very great accessibility in mouth



✓ 100% Dry air : the air is distributed by 4 conduits located around the hole of exit of water

✓ Air pressure allowing a fast drying in mouth as well as instruments or others

✓ Air / water mixture make inside of the cannula :

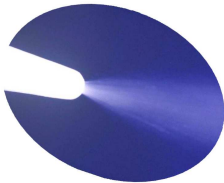

- a particularly homogeneous spray is obtained

- the power of the spray is increased by 30% compared to the original tip or with other disposable tips

✓ Adapters JMA is made out of anodized aluminium 20 µ: no deterioration of the appearance

Adapters availables : Luzzani, Siemens, Adec, Ritter top jet, Kavo, Faro, DCI, Colibri, Castellini, Chirana

